

General heuristic walkthrough to check on changes made since last evaluation.

Following rained then:

Description (last modified by wmoore) ([diff](#))

r: reporter o: owner

Must have priorities

General

- “Feed” label for home page
- Move the “Feed” button from the primary position ***DONE***
- User should first land on the ‘Data’ page after logging in ***DONE***
- Bell Icon for Activities

Replace this bell icon with a circular arrows icon to illustrate activities/active processes. See http://gretzky.openmicroscopy.org.uk/static/webgateway/img/activity_light_24.png)

Group

- As the ‘group’ dropdown is not always present, the main toolbar jumps around at the moment. Please can we show this button all the time even if inactive [Is this still an issue? If so, please give more details. Where do we need to show 'group' dropdown? (Will)]
- The dropdown menu spacing needs reduction
- The Tooltips interfere with use – please turn off [These actually have useful info: they say that you can click on the Group directly (don't **have** to choose a user). Maybe we could indicate this in another way? (Will)]
- Please order of the group members alphabetically **DONE**

History

- The calendar click through to image is not currently working

DONE?

- Filter appears in centre panel with the history – please remove this. [Is this not a useful feature? (Will)]
- The tabs of General, Acquisition and Preview should not be enabled initially. **DONE**

Web Admin

- There is no indication of whether you are currently on the ‘edit groups’ or ‘edit scientists’ tab. Please show a depression or colour contrast to illustrate this here. [Classes have been added, just needs css work (Will)]
- Web Admin Cog icon

Replace the cog icon in the web admin for editing scientists and groups with the pencil icon (this will achieve consistency across the application and will remove re-use of cogs on icons for different functionalities).

Search

- Only show Preview & Acquisition tabs when selecting an Image

DONE

- Only show search filter once we've done a search. **DONE**

Second Level Priorities

Right hand panel

- Accordion UI for Acquisition metadata should have some indication that it can be expanded **DONE**

General

- Thumbs look funny when loading with 'image' text and the whole page flowing as thumbs are loading. See screen-shots.
- The thumbnails/ dropdown looks like a button, but it actually works like a drop down menu. It should appear more like a drop-down menu - e.g. arrow icon on right hand side. **DONE**


- Profile icon defaults should include more diverse options. Alternatively just remove the default icon (Inform Insight)

- The numbers to the left of a project/dataset name are still at a low contrast level. Need to increase contrast and suggest using the brackets Insight uses i.e. [4] [You mean in the tree? I actually quite like it as it is (Will)]

To be checked:

Feed button position - much better visually with data as default home page

Bell icon - still there

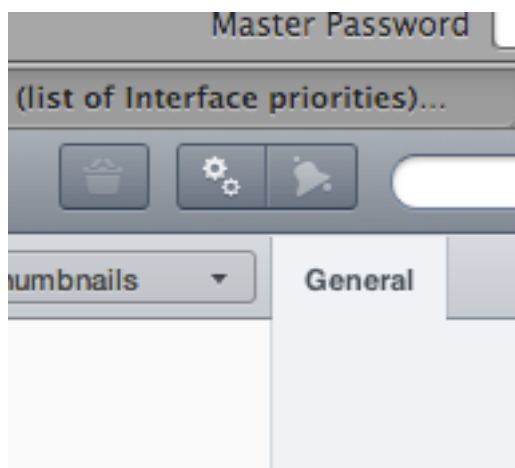

Group dropdown - sorted

Tooltips on group - better position but still a bit problematic

Order of group members - good much better

Calendar link - working

Filter - see below

Right hand tabs - good

Tabs group window - still and issue

Cog icon - replaced


Filter only after search - done

Aquisition expand


Thumbnails dropdown

Group:

Not sure why group coming up as "default"

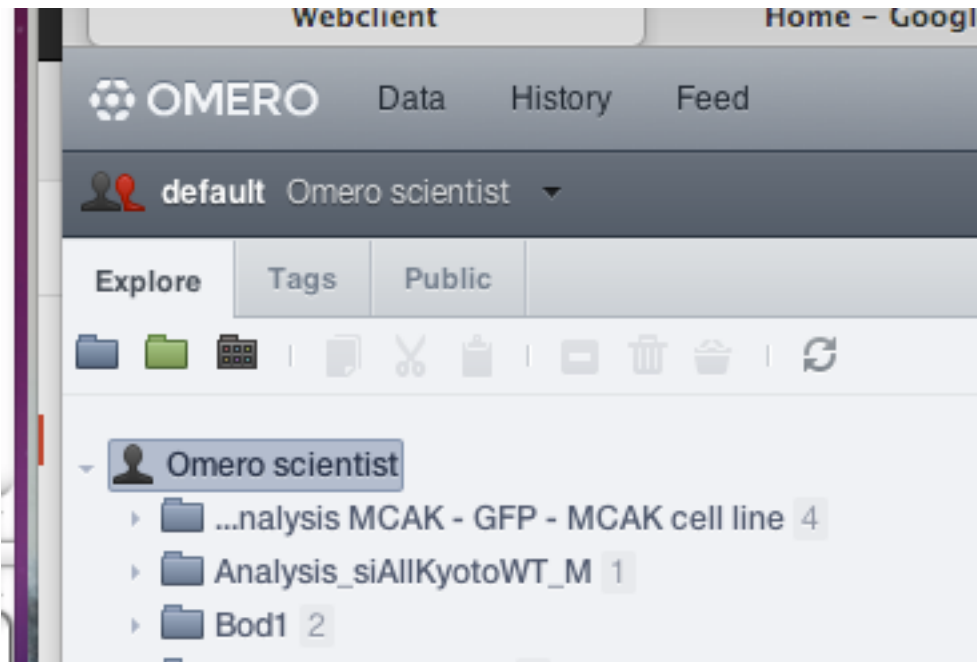

When I am a member of more than one group, I am still getting "default" Surely my default should be Swedlow

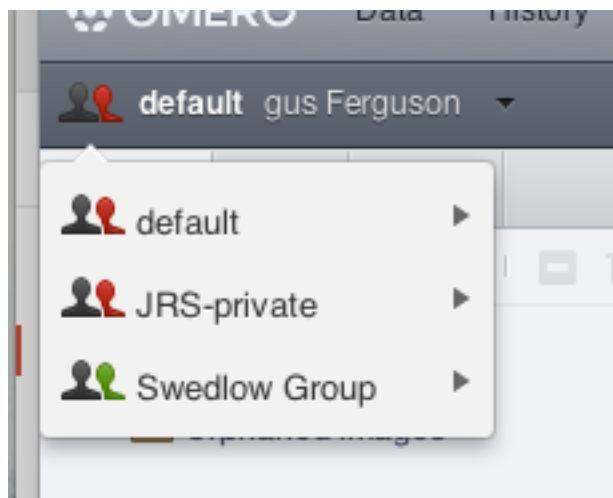



When I select omero scientist from the Swedlow group all his data is locked, and he was not a member of the Swedlow group when I was logged in as him. and had different data - are there 2 of him?



Bug Note:


In user account group “tab” remains available but is not active - should be greyed out if not active

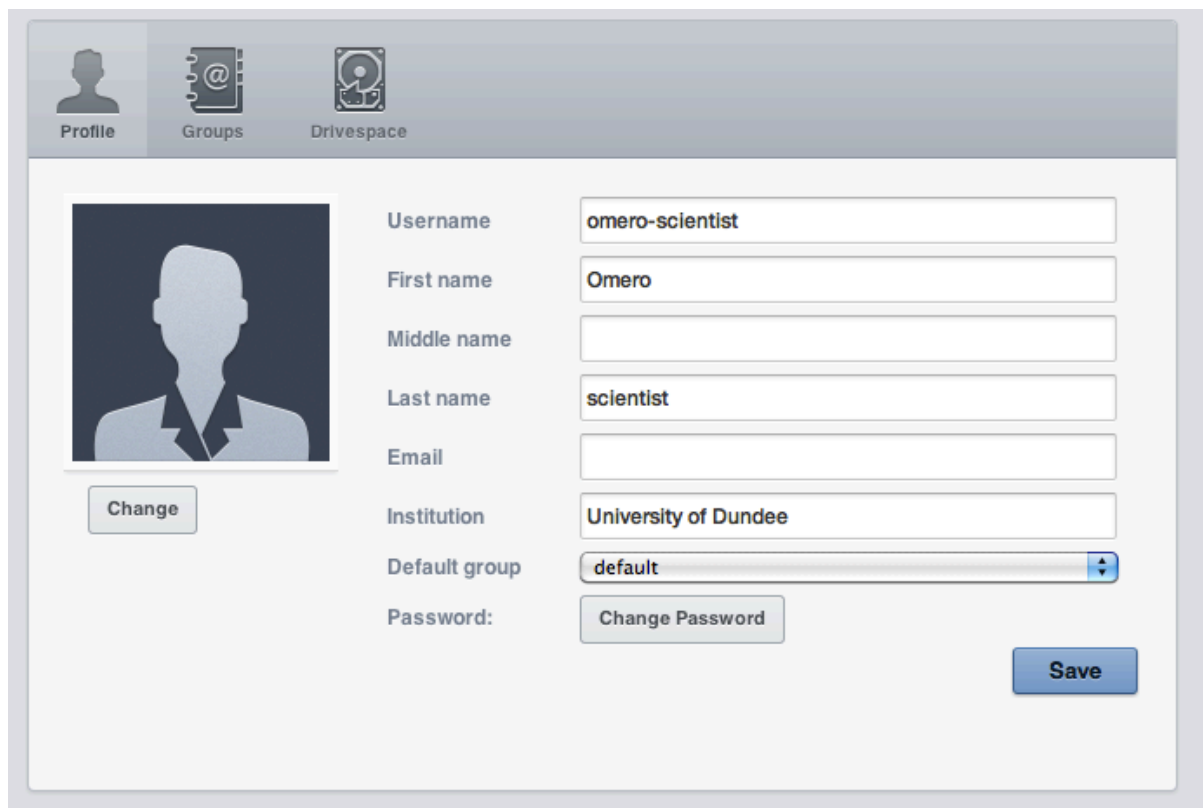

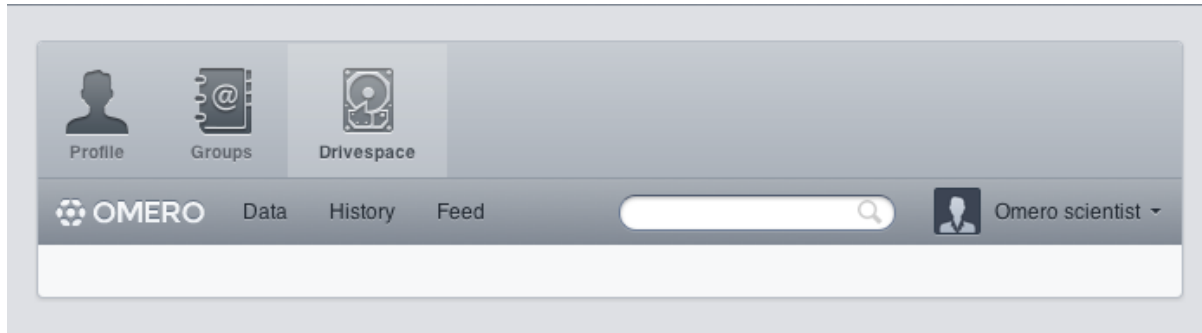

See separate bug document


When trying Drivespace I get my account but as dropdown with no functionality.

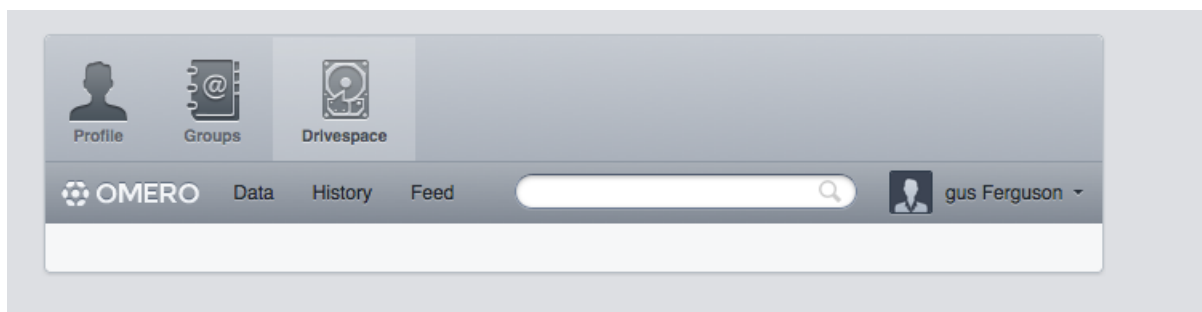


Bug Note:

History Tab

- get persistent loading icon by filter box when in History tab
- goes away when there is data in the centre pane

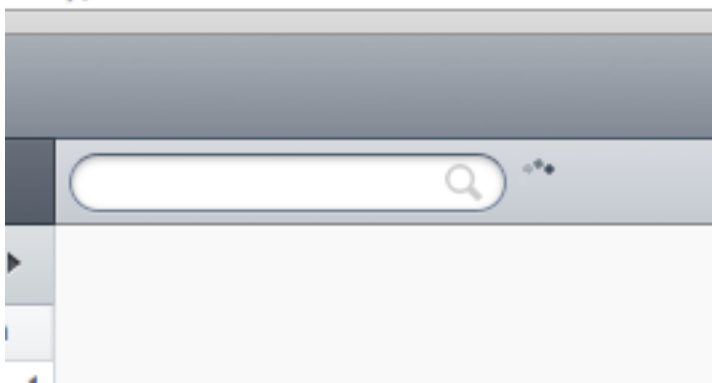


“Legend” box popping up at bottom of search results.
Pointless and out of context - suggest get rid of it all together

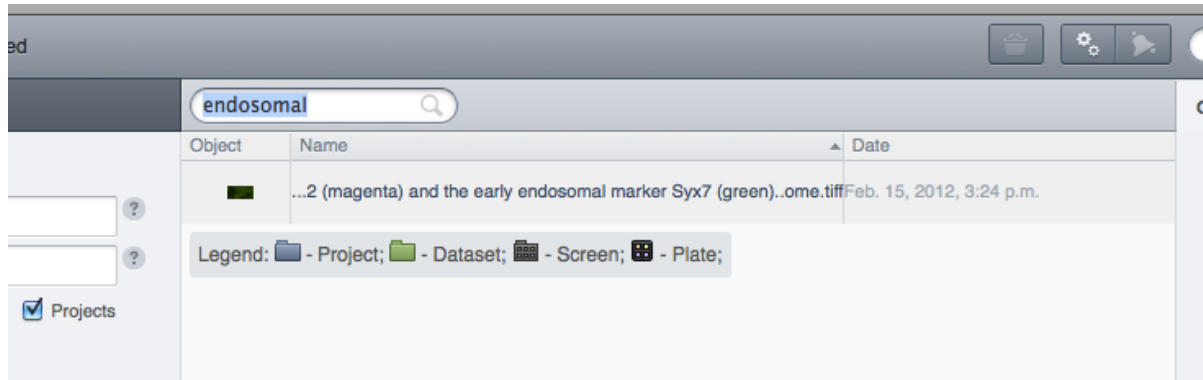


Bug Note:

Search Issues:


Transparent background to dates dropdown when clicked in Date box - needs solid colour BG

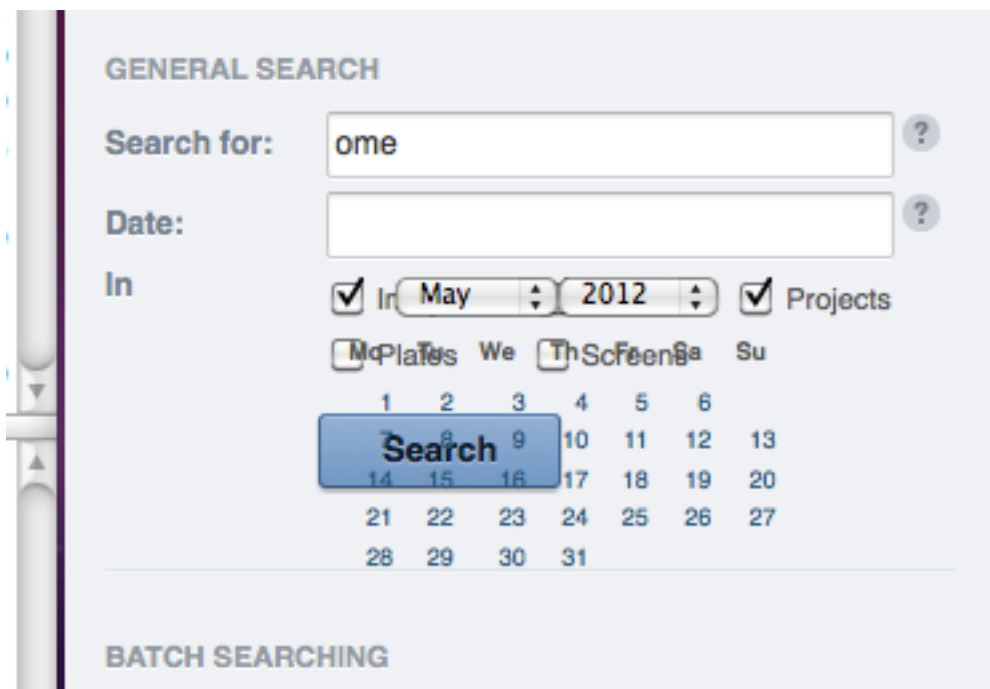


Date box label replaced by "In" when I pressed <cmd> to do screen capture - weird!

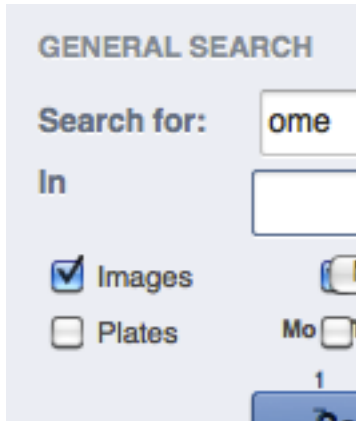


Search for “centrosomes” - failed to find any images when there is an image and a project with it in title

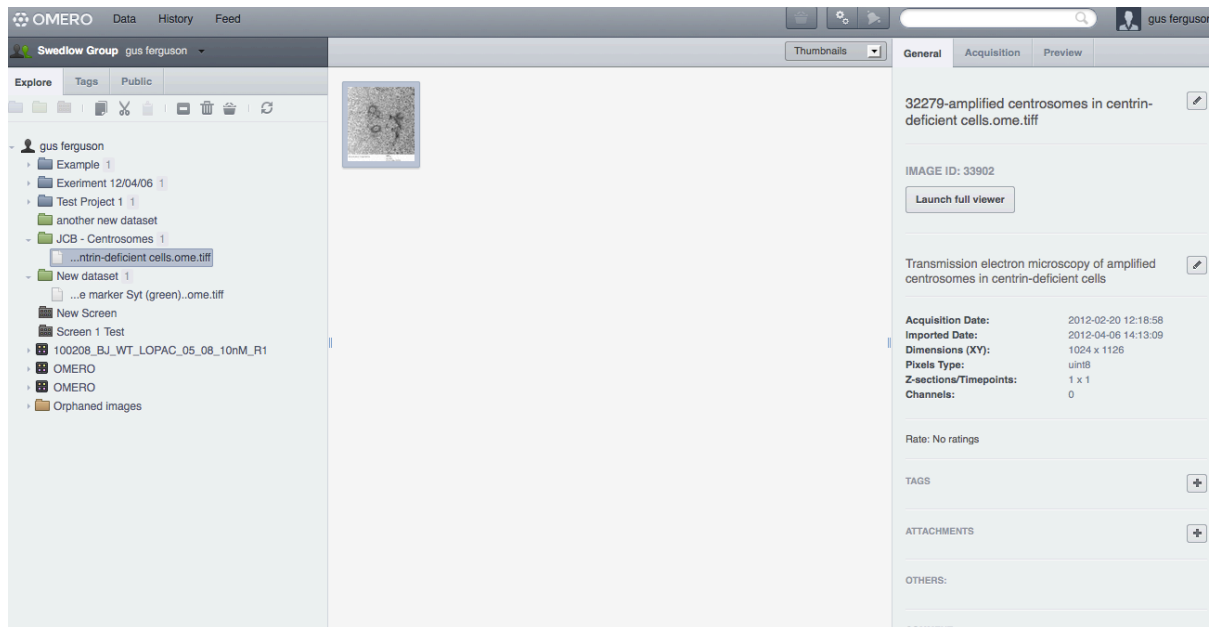

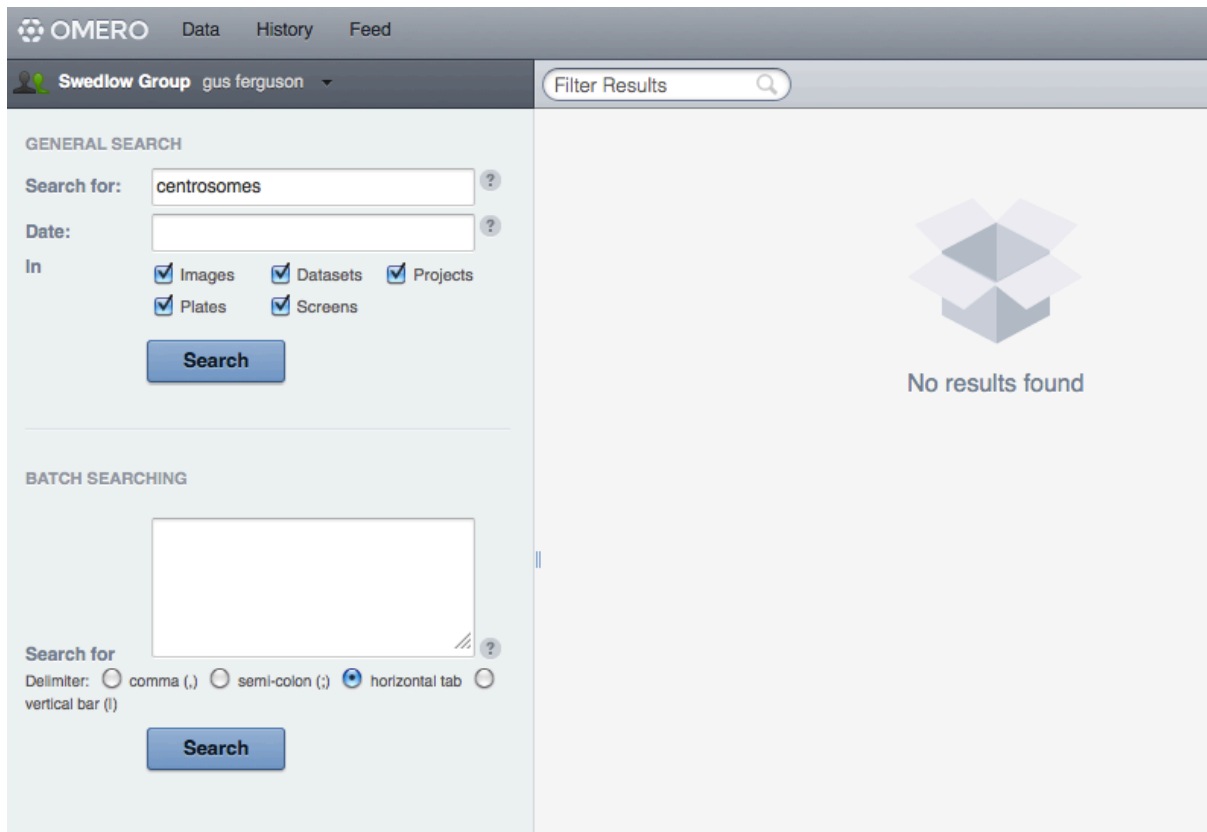

Search for "ome" - found images with no apparent "ome" in them

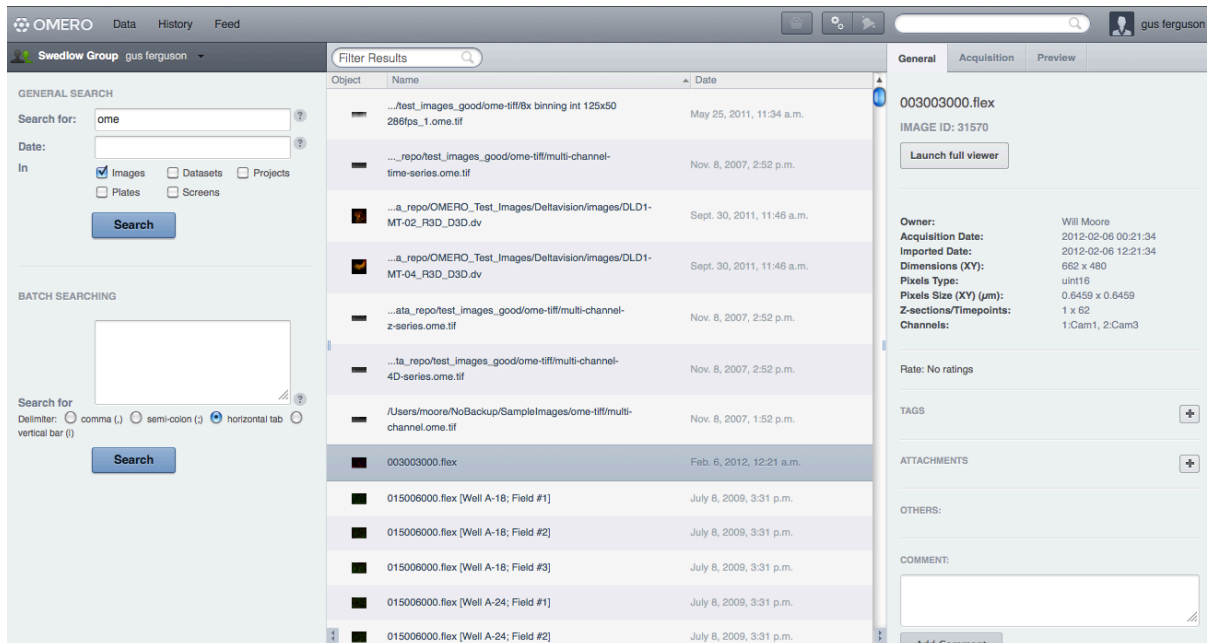

The screenshot shows the OMERO web interface. On the left, there is a search panel with a search box containing 'ome', a date field, and checkboxes for 'Images', 'Datasets', 'Projects', 'Plates', and 'Screens'. Below this is a 'BATCH SEARCHING' section with a larger search box and a 'Search' button. The main central pane displays a table of search results with columns for 'Object', 'Name', and 'Date'. The table lists various image files, including '015006000.flex [Well A-18; Field #1]'. On the right, a metadata panel for the selected image shows details such as 'Owner: Jason Swedlow', 'Acquisition Date: 2009-07-08 15:31:24', 'Imported Date: 2011-02-19 21:33:41', 'Dimensions (XY): 664 x 491', 'Pixels Type: uint16', and 'Z-sections/Timepoints: 1 x 1'. The metadata panel also includes sections for 'Rate', 'TAGS', 'ATTACHMENTS', 'OTHERS', and 'COMMENT'.

Inconsistent behaviour of filter box:

Did appear and stay in another search - see top screenshot

Filter box did not appear above results after another search - see screenshot above

No way of getting it back without repeating search

If I do search, get results in centre pane, then click on image, metadata for that image appears in right hand pane.

If I then do another search from left hand pane, I get new results in middle pane, but the metadata from the old image persists in right hand pane - needs refreshing with new search


Refreshes OK if new search done from toolbar search box.

Inconsistent search results:

Did search for "ome"

Got:

This a couple of times:

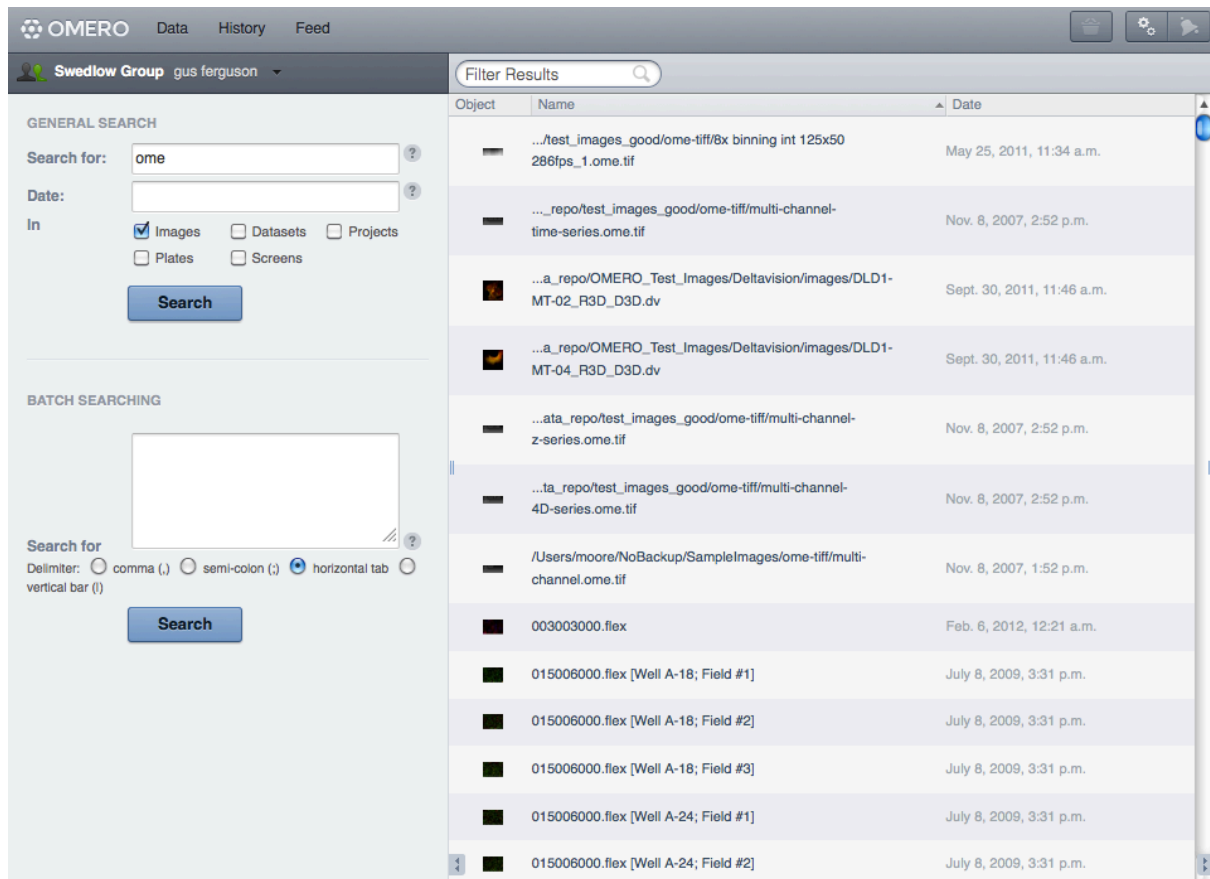

The screenshot shows the OMERO web interface. The top navigation bar includes 'Data', 'History', and 'Feed'. The user is logged in as 'gus ferguson' under the 'Swedlow Group'. The search results are displayed in a table with columns for 'Object', 'Name', and 'Date'. The search term 'ome' is entered in the 'GENERAL SEARCH' section. The results include various image and dataset files, such as '286fps_1.ome.tif' and '015006000.flex'.

Object	Name	Date
...	.../test_images_good/ome-tiff/8x binning int 125x50 286fps_1.ome.tif	May 25, 2011, 11:34 a.m.
...	..._repo/test_images_good/ome-tiff/multi-channel-time-series.ome.tif	Nov. 8, 2007, 2:52 p.m.
...	...a_repo/OMERO_Test_Images/Deltavision/images/DLD1-MT-02_R3D_D3D.dv	Sept. 30, 2011, 11:46 a.m.
...	...a_repo/OMERO_Test_Images/Deltavision/images/DLD1-MT-04_R3D_D3D.dv	Sept. 30, 2011, 11:46 a.m.
...	...ata_repo/test_images_good/ome-tiff/multi-channel-z-series.ome.tif	Nov. 8, 2007, 2:52 p.m.
...	...ta_repo/test_images_good/ome-tiff/multi-channel-4D-series.ome.tif	Nov. 8, 2007, 2:52 p.m.
...	/Users/moore/NoBackup/SampleImages/ome-tiff/multi-channel.ome.tif	Nov. 8, 2007, 1:52 p.m.
...	003003000.flex	Feb. 6, 2012, 12:21 a.m.
...	015006000.flex [Well A-18; Field #1]	July 8, 2009, 3:31 p.m.
...	015006000.flex [Well A-18; Field #2]	July 8, 2009, 3:31 p.m.
...	015006000.flex [Well A-18; Field #3]	July 8, 2009, 3:31 p.m.
...	015006000.flex [Well A-24; Field #1]	July 8, 2009, 3:31 p.m.
...	015006000.flex [Well A-24; Field #2]	July 8, 2009, 3:31 p.m.

then suddenly got this:

No reason for different results apparent

OMERO Data History Feed

Swedlow Group gus ferguson

Filter Results

GENERAL SEARCH

Search for:

Date:

In Images Datasets Projects
 Plates Screens

Search

BATCH SEARCHING

Search for

Delimiter: comma (,) semi-colon (;) horizontal tab vertical bar (|)

Search

Object	Name	Date
	.../test_images_good/ome-tiff/8x binning int 125x50 286fps_1.ome.tif	May 25, 2011, 11:34 a.m.
	..._repo/test_images_good/ome-tiff/multi-channel-time-series.ome.tif	Nov. 8, 2007, 2:52 p.m.
	...a_repo/OMERO_Test_Images/Deltavision/images/DLD1-MT-02_R3D_D3D.dv	Sept. 30, 2011, 11:46 a.m.
	...a_repo/OMERO_Test_Images/Deltavision/images/DLD1-MT-04_R3D_D3D.dv	Sept. 30, 2011, 11:46 a.m.
	...ata_repo/test_images_good/ome-tiff/multi-channel-z-series.ome.tif	Nov. 8, 2007, 2:52 p.m.
	...ta_repo/test_images_good/ome-tiff/multi-channel-4D-series.ome.tif	Nov. 8, 2007, 2:52 p.m.
	/Users/moore/NoBackup/SampleImages/ome-tiff/multi-channel.ome.tif	Nov. 8, 2007, 1:52 p.m.
	COLUMBUS.jpg	Aug. 21, 2011, 3:05 p.m.
	Face.jpg	Aug. 29, 2011, 11:10 a.m.
	import.jpg	Aug. 21, 2011, 3:36 p.m.
	JCB DV.jpg	Aug. 21, 2011, 4:13 p.m.
	JCB_DV.png	Aug. 21, 2011, 4:16 p.m.
	multi-channel	Nov. 8, 2007, 1:52 p.m.

BUG: - Tried to select another user under Groups menu:

Using Safari.

This had worked fine earlier.

Same happened with trying to select any other users.

Failed to replicate in Firefox with Firebug on.

Go back to: [Webadmin](#), [Webclient](#) or [Logout](#)

Server Error. (500)

There's been an error. Please give us feedback to help to improve the OMERO applications. Thanks for your patience.

Your email

Comment

Traceback (most recent call last):

```

File
"/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/core/handlers/base.py", line 111, in get_response
    response = callback(request,
*callback_args, **callback_kwargs)

```

Error

{{

Traceback (most recent call last):

```

File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/core/handlers/base.py", line 111, in get_response
    response = callback(request, *callback_args, **callback_kwargs)

```

```

File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/omeroweb/decorators.py", line 313, in wrapped
    retval = f(request, *args, **kwargs)

```

```

File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/omeroweb/decorators.py", line 347, in wrapper
    context = f(request, *args, **kwargs)

```

```

File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/omeroweb/webclient/views.py", line 368, in load_template
    if not form_users.is_valid():

```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/forms/forms.py", line 121, in is_valid
    return self.is_bound and not bool(self.errors)
```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/forms/forms.py", line 112, in _get_errors
    self.full_clean()
```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/forms/forms.py", line 267, in full_clean
    self._clean_fields()
```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/forms/forms.py", line 284, in _clean_fields
    value = field.clean(value)
```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/django/forms/fields.py", line 163, in clean
    value = self.to_python(value)
```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/omeroweb/webadmin/custom_forms.py", line 331, in to_python
    if checkValue(qq, value):
```

```
File "/home/omero/slave/workspace/OMERO-merge-green/src/dist/lib/python/omeroweb/webadmin/custom_forms.py", line 323, in checkValue
    if long(value) == q.id:
```

ValueError: invalid literal for long() with base 10: '4?Acquisition=1101'

<WSGIRequest

GET: <QueryDict: {u'experimenter': [u'4?Acquisition=1101']}>,

POST: <QueryDict: {}>,

COOKIES: {'sessionid': '91dcc02e17983a68ce4aa9ce9e535c5c'},

META: {'DOCUMENT_ROOT': '/var/www',

'GATEWAY_INTERFACE': 'CGI/1.1',

'HTTP_ACCEPT': 'text/html,application/xhtml+xml,application/xml;q=0.9,*/*;q=0.8',

'HTTP_ACCEPT_ENCODING': 'gzip, deflate',

'HTTP_ACCEPT_LANGUAGE': 'en-us',

'HTTP_CONNECTION': 'keep-alive',

'HTTP_COOKIE': 'sessionid=91dcc02e17983a68ce4aa9ce9e535c5c',

'HTTP_HOST': 'gretzky.openmicroscopy.org.uk',

'HTTP_REFERER': 'http://gretzky.openmicroscopy.org.uk/webclient/userdata/?experimenter=1553',


```
'HTTP_USER_AGENT': 'Mozilla/5.0 (Macintosh; Intel Mac OS X 10_6_8)
AppleWebKit/534.55.3 (KHTML, like Gecko) Version/5.1.5 Safari/534.55.3',
'PATH': '/usr/local/bin:/usr/bin:/bin',
'PATH_INFO': u'/webclient/userdata/',
'PATH_TRANSLATED': '/home/omero/OMERO-CURRENT/var/omero.fcgi/
webclient/userdata/',
'QUERY_STRING': 'experimenter=4?Acquisition=1101',
'REMOTE_ADDR': '10.34.1.133',
'REMOTE_PORT': '59044',
'REQUEST_METHOD': 'GET',
'REQUEST_URI': '/webclient/userdata/?experimenter=4?Acquisition=1101',
'SCRIPT_FILENAME': '/home/omero/OMERO-CURRENT/var/omero.fcgi',
'SCRIPT_NAME': u'',
'SERVER_ADDR': '134.36.65.227',
'SERVER_ADMIN': 'webmaster@localhost',
'SERVER_NAME': 'gretzky.openmicroscopy.org.uk',
'SERVER_PORT': '80',
'SERVER_PROTOCOL': 'HTTP/1.1',
'SERVER_SIGNATURE': '<address>Apache/2.2.14 (Ubuntu) Server at
gretzky.openmicroscopy.org.uk Port 80</address> \n',
'SERVER_SOFTWARE': 'Apache/2.2.14 (Ubuntu)',
'wsgi.errors': <flup.server.fcgi_base.TeeOutputStream object at
0x2f67d10>,
'wsgi.input': <flup.server.fcgi_base.InputStream object at 0x2c78ad0>,
'wsgi.multiprocess': True,
'wsgi.multithread': False,
'wsgi.run_once': False,
'wsgi.url_scheme': 'http',
'wsgi.version': (1, 0)}>
```

```
}}}
```

BUG: Uncaught Exception when Clicking on Groups Tag in “Account Settings”

Console:

uncaught exception: jQuery UI Tabs: Mismatching fragment identifier.

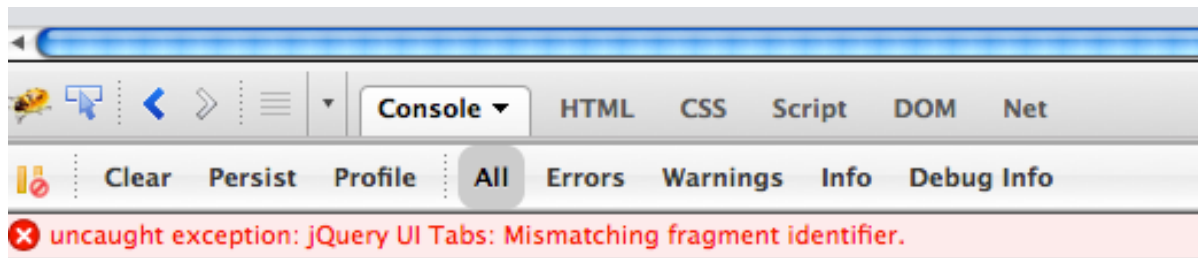

```
{{{
```

```
<!DOCTYPE html PUBLIC "-//W3C//DTD XHTML 1.0 Transitional//EN" "http://
www.w3.org/TR/xhtml1/DTD/xhtml1-transitional.dtd">
<html xmlns="http://www.w3.org/1999/xhtml">
<head>
<link type="image/x-icon" href="/favicon.ico" rel="shortcut icon">
<link type="image/x-icon" href="/favicon.ico" rel="icon">
<link type="text/css" href="/static/webgateway/css/reset.css" rel="stylesheet">
<link type="text/css" href="/static/webgateway/css/ome.body.css" rel="stylesheet">
<link type="text/css" href="/static/3rdparty/jquery-ui-1.8.19/themes/base/
jquery.ui.base.css" rel="stylesheet">
<link media="screen" type="text/css" href="/static/webclient/css/dusty.css" rel="stylesheet">
<link type="text/css" href="/static/webgateway/css/ome.header.css" rel="stylesheet">
<link media="screen" type="text/css" href="/static/webclient/css/
layout.css" rel="stylesheet">
<link media="screen" type="text/css" href="/static/webadmin/css/
layout.css" rel="stylesheet">
<link media="screen" type="text/css" href="/static/webgateway/css/
ome.table.css" rel="stylesheet">
<link media="screen" type="text/css" href="/static/3rdparty/jquery.tablesorter/
jquery.tablesorter.css" rel="stylesheet">
<title> My account </title>
<script src="/static/3rdparty/jquery-1.7.2.js" type="text/javascript">
<script>
<script src="/static/3rdparty/jquery-ui-1.8.19/jquery-ui-1.8.19.custom.min.js" type="text/
javascript">
<script src="/static/webgateway/js/ome.popup.js" type="text/javascript">
<script>
<script>
<script src="/static/3rdparty/jquery.tablesorter/jquery.tablesorter.js" type="text/javascript">
<script src="/static/3rdparty/jquery.quicksearch.js" type="text/javascript">
<script src="/static/webclient/javascript/jquery.form.js" type="text/javascript">
<script type="text/javascript">
<script>
<script>
<link type="image/x-icon" href="/favicon.ico" rel="shortcut icon">
</head>
<body>
<div id="middle_header">
```

```

<div id="content">
<div class="ui-dialog ui-widget ui-widget-content ui-corner-all reddit" style="display: none;
z-index: 1000; outline: 0px none;" tabindex="-1" role="dialog" aria-labelledby="ui-dialog-
title-password_change_dialog">
</body>
</html>

```

```
}}}
```

BUG: Syntax error

Opened in Image Viewer using the "Launch Full Viewer" button in the General tab of the right hand pane.

Syntax error caught in Firebug output.

Reproducible with other images.

```
{{{
```

```

GET http://gretzky.openmicroscopy.org.uk/webgateway/... 200 159ms
205483483888307329_1335966292576&_1335966293006 OK ms
jquery-1.7.2.js
(line 8240)

```

Syntax error, unrecognized expression: [value=0.9]

```

    throw new Error( "Syntax error, unrecognized
expression: " + msg );
jquery-1.7.2.js (line 4267)

```

```
}}}
```



003003000.flex

gretzky.openmicroscopy.org.uk/webgateway/img_detail/31570/

Viewing Options

- Normal
- Max Intensity
- Split Channel
- Quality: Normal
- Zoom: 66
- Line Plot
- Rendering Details**
- Channels - Edit
- Color
- Current Image**
- Z: ?? | T: 1/1
- [Image Information](#)
- [Image Link](#)
- ROIs: [Show](#) | [Hide](#)

Z-sections

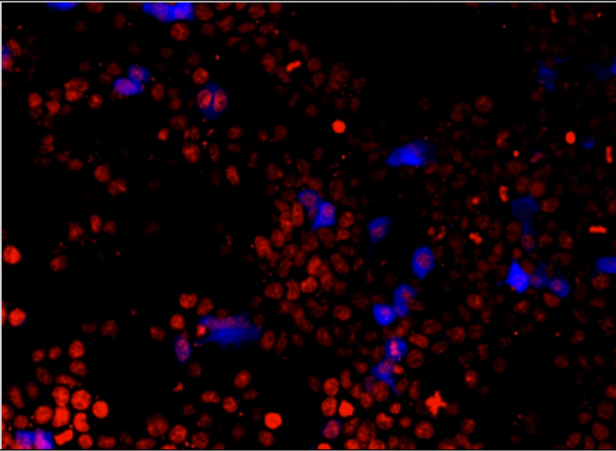

Timepoints

© 2007-2011 Glencoe Software Inc. All rights reserved.

Console HTML CSS Script DOM Net

Clear Persist Profile All Errors Warnings Info Debug Info

```
GET http://gretzky.openmicroscopy.org.uk/webgateway/...205483483888307329_1335966292576&_=133596629300 jquery-1.7.2.js (line 8240)  
Syntax error, unrecognized expression: [value=0.9] jquery-1.7.2.js (line 4267)  
throw new Error( "Syntax error, unrecognized expression: " + msg );
```

>>>

1P